

## Tracoe Twist Extract – REF 888-306-P

Tracheostomy tube with minimally traumatic inserter, low pressure cuff, subglottic suction and fenestration

### Product Description

The Tracoe Twist tracheostomy tubes are used as ventilation and therapy tubes with inner cannulas in hospitals, rehabilitation, and community settings. A special feature is the curved neck flange that moves in both the vertical and horizontal axes, following the movements of the patient's head. The thin tube walls of inner and outer cannulas allow for a good inside - outside diameter ratio. Tracoe Twist and Twist Plus tracheostomy tubes are made of polyurethane (PU). The low-pressure cuff is made of medical grade PVC.

### Scope of Supply

- Inner cannula, fenestrated with 15 mm connector
- 2 inner cannulas with 15 mm connector
- Occlusion cap
- 2 connectors for suction devices
- Reinsertion tube
- Lubricating gel
- Perforated obturator
- Neck strap




General Information		Material Specifications	
Product line	Twist Extract	Material Outer cannula	Polyurethane
Legal manufacturer	Tracoe Medical GmbH Reichelsheimer Str. 1/3 55268 Nieder-Olm, Germany	Material Neck Flange	Polycarbonate
Country of origin	Germany	Material Cuff	Polyvinylchloride
Colour Coding	REF 888-306-P	Biocompatibility	EN ISO 10993-1
Handling			
Shelf life	5 years, from date of production	Storage	
Maximum Useful Life	29 days, from first use		
Maximum Useful Life Insertion System	single-use disposable product		
Packaging		Certification	
Packaging unit	1 unit	Product classification acc. MDR	IIb
Weight	See Table	CE Marking	
Dimensions (W x H x L)	271 x 47 x 113	GMDN	35404
Material outer packaging:	Cardboard box	Quality system certified to	EN ISO 13485
Material sterile packaging:	Sterile bag with medical paper and transparent heatsealable PET/PP foil	FDA registration	n. A.

## Product Range

Size	Order Number	ID-IC	OD-OC bc	OD-OC bf	TL	OB	BA	CD	Gross Weight	Global Trade Item Number
		mm	mm	mm	mm	mm	°	mm	g	
07	REF 888-306-07-P	7,0	10,4	12,0	74,0	80,5	95	26,0	236	4035324019120
08	REF 888-306-08-P	8,0	11,4	12,7	76,0	83,0	90	28,0	239	4035324019137
09	REF 888-306-09-P	9,0	12,5	14,2	78,0	88,5	90	30,0	242	4035324019144

## Legend

- ID-IC** Inside diameter at bottom of inner cannula
- OD-OC bc** Outside diameter at bottom of outer cannula
- OD-OC bf** Outside diameter of outer cannula behind the flange
- TL** Length along centre line from start of neck flange to bottom of tube
- OB** Length along outer bend from start of neck flange to bottom of tube
- BA** Bending angle
- CD** Cuff diameter

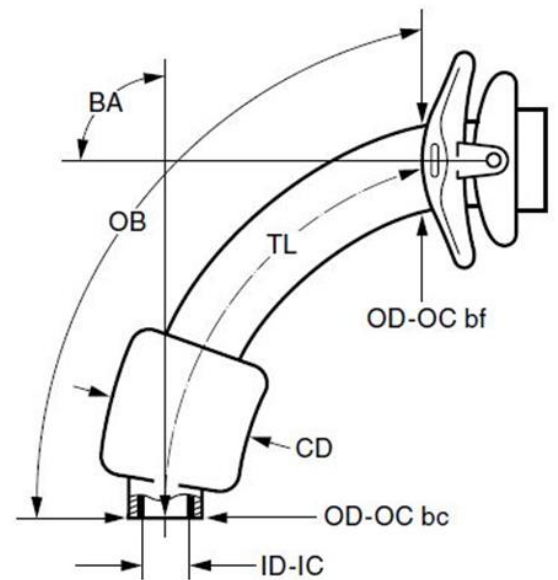

Breathing-Speaking-Living

Tracoe Medical GmbH  
Reichelsheimer Str. 1 / 3  
55268 Nieder-Olm  
Germany

Atos Medical AB  
Hyllie Boulevard 17  
215 32 Malmö  
Sweden  
[www.atosmedical.com](http://www.atosmedical.com)