

Tracoe Twist Plus – REF 314

Tracheostomy tube with double fenestration at inner and outer bend

Product Description

The Tracoe Twist Plus tracheostomy tubes are used as ventilation and therapy tubes with inner cannulas in hospitals, rehabilitation, and community settings. A special feature is the curved neck flange that moves in both the vertical and horizontal axes, following the movements of the patient's head. The thin tube walls of inner and outer cannulas allow for a good inside - outside diameter ratio. Tracoe Twist and Twist Plus tracheostomy tubes are made of polyurethane (PU). The low-pressure cuff is made of medical grade PVC.

Scope of Supply

- 1 Inner cannulas with 15 mm connector
- 1 Inner cannula, double fenestrated with 15 mm connector
- Occlusion cap
- Perforated obturator
- Neck strap



General Information	
Product line	Twist Plus
Legal manufacturer	Tracoe Medical GmbH Reichelsheimer Str. 1/3 55268 Nieder-Olm, Germany
Country of origin	Germany
Colour Coding	REF 314 ++
Handling	
Shelf life	5 years, from date of production
Maximum Useful Life	5 weeks, from first use

Packaging	
Packaging unit	1 unit
Weight	See Table
Dimensions (W x H x L)	210 x 113 x 50 mm
Material outer packaging:	Cardboard box
Material sterile packaging:	Sterile bag with medical paper and transparent heatsealable PET/PP foil

Material Specifications	
Material Outer cannula	Polyurethane
Material Neck Flange	Polycarbonate
Material Cuff	n. A.
Biocompatibility	EN ISO 10993-1
	

Storage

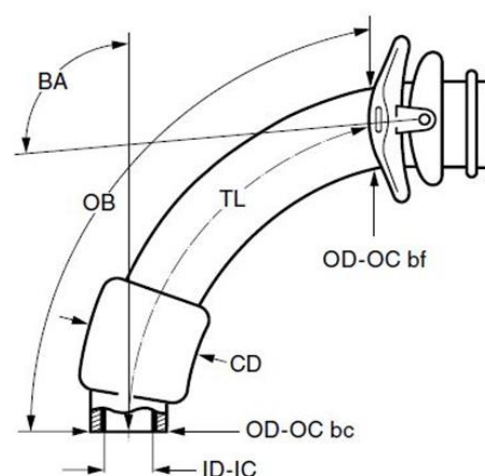

Certification	
Product classification acc. MDR	IIb
CE Marking	
GMDN	35404
Quality system certified to	EN ISO 13485
FDA registration	Clearance Number K961449

Product Range

Size	Order Number	ID-IC	OD-OC bc	OD-OC bf	TL	OB	BA	Gross Weight	Global Tarde Item Number
		mm	mm	mm	mm	mm	°	G	
07	REF 314-07	7,0	9,8	10,1	85	91	100	110	4035324018154
08	REF 314-08	8,0	10,8	11,1	88	95	100	110	4035324018161
09	REF 314-09	9,0	11,8	12,1	90	99	100	112	4035324018178
10	REF 314-10	10,0	12,8	13,1	92	102	100	113	4035324018185

Legend

ID-IC	Inside diameter at bottom of inner cannula
OD-OC bc	Outside diameter at bottom of outer cannula
OD-OC bf	Outside diameter of outer cannula behind the flange
TL	Length along centre line from start of neck flange to bottom of tube
OB	Length along outer bend from start of neck flange to bottom of tube
BA	Bending angle
CD	Cuff diameter



Tracoe Medical GmbH
Reichelsheimer Str. 1 / 3
55268 Nieder-Olm
Germany



Atos Medical AB
Hyllie Boulevard 17
215 32 Malmö
Sweden
www.atosmedical.com