




Product Information

TRACOE Connection Tube with check valve

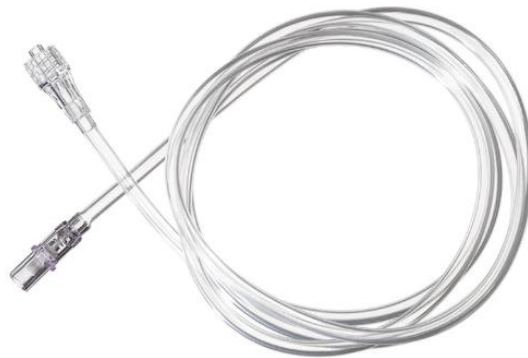

Product description:

The **TRACOE Connection Tube with check valve** is a non-sterile, 1 m long Luer-to-Luer tube designed to connect a cuff pressure monitor to a tracheostomy or endotracheal tube cuff system. One end of the TRACOE Connection Tube connects to the pilot balloon valve (Luer-compatible). The other end connects to a monitor or inflation device (Luer-compatible).

The device can be used on adults, pediatrics, and neonates for a maximum of 72 hours in total.

Manufacturer: PROMEPLA S.A.M.
9 Avenue Albert II
98000 Monaco

Classification: Class I (rule 1)
(EU) 2017/745

Intended Use: The **TRACOE Connection Tube with check valve** is designed for the pneumatic connection between the *TRACOE Cuff Pressure Monitor Sensitive/TRACOE Smart Cuff Manager* and the valve of the cuff inflation line/pilot balloon of a tracheostomy tube or an



CE Mark:	endotracheal tube. Yes, the device is CE marked.
EMDN code :	A03020101 - Low pressure extension lines
Sterilization:	Non-sterile
Raw material:	PVC (<i>Polyvinyl Chloride</i>), PBT (<i>Polybutylene Terephthalate</i>), NBR (<i>Nitrile Butadiene Rubber</i>), Stainless steel, ABS (<i>Acrylonitrile Butadiene Styrene</i>)
Latex information:	Not manufactured with natural rubber latex.
Biological origin:	The device is not manufactured with materials derived from human or animal source.
Handling and storage:	Store the product dry and away from sunlight.
Waste handling and disposal:	Waste handling and disposal should be carried out in agreement with medical practice and applicable national laws and legislations.
Hazardous components:	None
Expiration date:	5 years after manufacturing.
Packaging:	10 pcs in a box, each individually packed in a peel pouch (PET/PP/PAPER)
Length:	1 meter
Distributed by :	Coloplast A/S Holtedam 1 3050 Humlebaek Denmark

Devices under Basic UDI-DI: 037005129ACC06N2

REF	Name	UDI-DI
ODG0702SC167	TRACOE Connection Tube with check valve	03700512906995

Atos Medical AB Compatible products:

Name	Ref No.
<i>TRACOE Cuff Pressure Monitor Sensitive</i>	REF 721
<i>TRACOE Smart Cuff Manager</i>	REF 730-5